


Flowchart

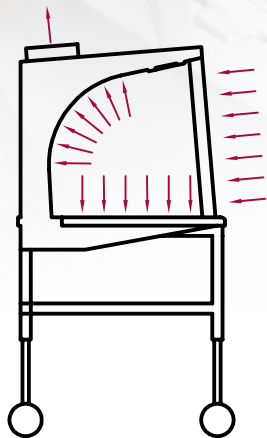

Options:

- Anesthesia induction box
- Weighing section
- Power and data output
- Pressure and flow alarms

Ventilated table for operator protection

Ventilated table for operator protection and allergen control. Ideal for cage change, procedures, dosing, anesthesia...

The air is sucked in at the front opening and down through the perforated table top and the rounded rear wall.

Electrically height adjustable. Integrated and dimmable LED lighting. G4 prefilter in rear wall.

For connection to room exhaust. Mobile on wheels or on fixed legs.

Designed in cooperation with JRV A/S – a lifetime of experience within air flow solutions



Specifications:

- **Top:** Stainless steel or Corian (white) and transparent plastic side panels in PTEG
- **Sink:** Stainless steel
- **Frame:** 40/40 coated square tubing RAL 9010 white
- **Controls:** Operator control panel on front, dimmable LED light, visual pressure alarm
- **Ext. Dim:** W X D X H: 1800 / 1600 / 1200 x 780 x 720 – 1100mm
- **Int. dim:** W X D X H: 1740 / 1540 / 1140 x 555 x 518mm
- **Power:** 230 volt, 50 Hz
- **Exhaust connection to room ventilation:** Ø250mm

Available in 1200, 1600 and 1800 wide models. Required air volume depends on the size.